

# Water Baths

## Digital Series



The SHEL LAB high-performance water baths are accurate, easy to use, safe and durable. The water bath design incorporates a drip free cover holster and pocket handles so users can easily transport the water bath.

A microprocessor achieves precise temperature control regardless of how the unit is loaded. Calibration is performed with the convenient, front panel touch pad, which is essential for facilities concerned with IQ/OQ/PQ validation.

### FEATURES

- Easy-to-Clean Stainless Steel Tank
- Pocket Handles for Easy Lifting
- Recessed Heating Element Prevents "Burnout"
- Unique SHEL LAB Design Eliminates "Hot Spots"
- Thermometer Clip Included

SHEL LAB was the first to introduce the non-contact recessed heating element design to the analytical research market. This curtails element burnout and eliminates tank hot spots that are chronic challenges for other water baths.

The SHEL LAB digital water baths remain an industry standard for precise control and quality of design. The SHEL LAB constant temperature digital water baths are perfect for conducting a host of applications including: bacteriological examinations, food processing/QC procedures and microbiology assays just to name a few.

*SHEL LAB digital water baths are durable and dependable.*

### Water Bath Applications

- Sample Thawing
- Bacteriological Examinations
- Warming Reagents
- Coliform Determinations
- Microbiological Assays



## Precise Temperature Control -

- Over Temperature Protection
- Temperature Uniformity +/- 0.2% at 37°C
- Temperature Range Ambient + 5°C to 80°C
- Coated tanks available for applications up to 80°C




General Purpose Water Baths	110-120V 220-230V Chambers	Model Number		
		SWB7 SWB7-2 Single (1)	SWB15 SWB15-2 Single (1)	SWB23 SWB23-2 Single (1)
Exterior Dimensions (wxdxh)	Inches	15.0 x 12.5 x 12.0	15.0 x 19.0 x 12.0	14.8 x 25.0 x 12.0
	mm	381 x 318 x 305	381 x 483 x 305	376 x 635 x 305
Bath Tank Dimensions (wxdxh)	Inches	11.8 x 5.8 x 6.0	11.8 x 13.0 x 6.0	11.5 x 19.5 x 6.0
	mm	300 x 174 x 152	300 x 330 x 152	292 x 495 x 152
Tank Capacity	L	6	14	23
Stainless Steel Gable Cover(s)	Yes/No	Yes - One Included	Yes - One Included	Yes - One Included

\*Tank is excluded from the 2 year warranty. All specifications are determined by using average values on standard equipment at an ambient temperature of 25°C (77°F) and line voltages within +/-10% of unit type (115V/230V). Temperature specifications follow DIN 12880 methodology. We reserve the right to change specifications at any time.



0740514 3/14

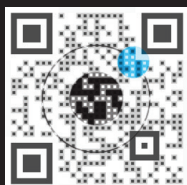

www.shellab.com  
1-800-322-4897